

HPLC BENCH BARCELONA

Keeping your equipment
in good hands.


.....



Keeping your equipment
in good hands.



Scan the QR code
to see this product
on our website or
[click here](#)

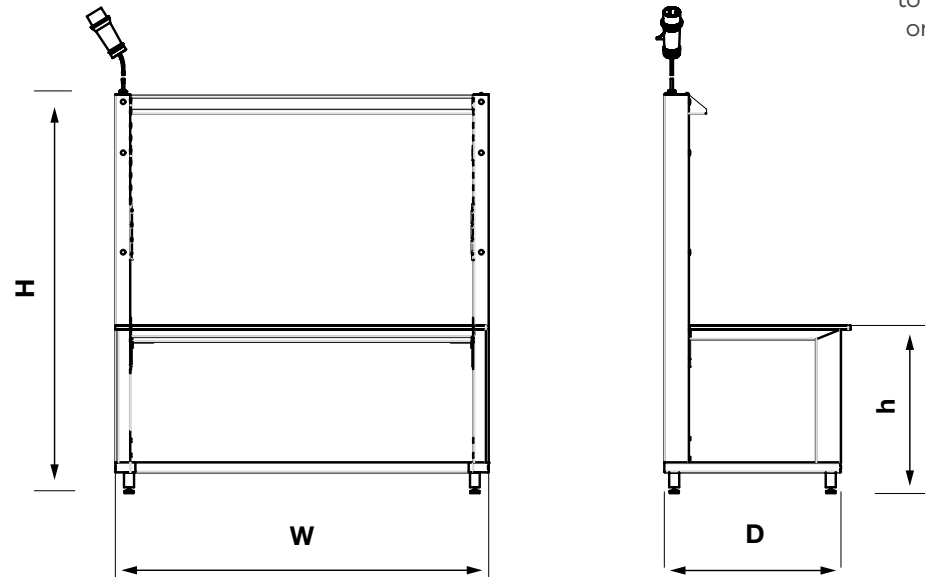

W (width, mm)	D (depth, mm)	h (worktop height, mm)	H (column height, mm)
1500			
1200	750	750	1570
900	900	650	1800
750			

MAX Load 350 Kg

Power supply side can be switched left to right and viceversa. Check **accessories pages** to see how shelves can be added, removed and changed position.

- **Underbench shelf** available (highly advised for waste collection).
- **Gas taps** can be attached to column.
- Some **other options** available.
- In addition, **special accessory attachments** can be design and manufactured on demand.

Main frame fully welded steel, epoxy coated 9010 (check other options) with M5 nuts below worktop for accesories and attachments.

Columns manufactured in steel tube bolted to main frame, epoxy coated 9010 (check other options), with **M8 nuts at internal faces** for accesories. **Nuts in front faces** as an option.

Worktop material on demand. Advised HPL 16 mm high chemical resistance.

Feet, articulated with rubber base, interchangeable with combi castor-foot.

Data setup on demand, advised to be magnet clip-routed.

Plug + Play interchangeable set-up.

Typical connection to ceiling options:

- CE connector plug
- Schuko plug
- UK plug
- Thorsman plug

M8 capped nuts for accesories & attachments, at front face of columns.


(x3)

Bench can be supplied with a circuit breaker set mounted below sockets at left column (custom setup per project)

Leveling feet and Combi castor-foot are interchangeable



Leveling feet

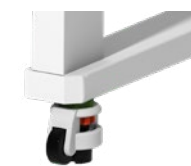


Combi castor-foot

Worktop is delivered with drills at front and back to allow waste piping routing to the bottom.

M8 Capped nuts for Shelving and Accesories, 70x50mm Pattern.

Quantity depends on column height.

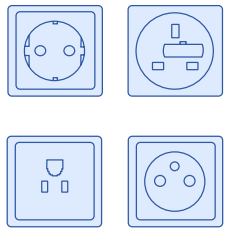

Most standard outlets available. 8 outlets per workstation: 3 outlets at each side above worktop. 1 outlet at each side below worktop.

Special short shelf with retaining bar allows for tall equipment to be installed while keeping items tight and secure.



Some examples of **worktop materials**.

(12x) M5 nuts provided below worktop, in pairs distanced 95 mm between axis (distance between pairs depends on bench size).