

# PRIMARY BENCH

*A clean, adaptable, and efficient  
solution for everyday lab work.*

.....




# Primary Bench 4P

A clean, adaptable, and efficient solution for everyday lab work.



Scan the QR code to see this product on our website or [click here](#)

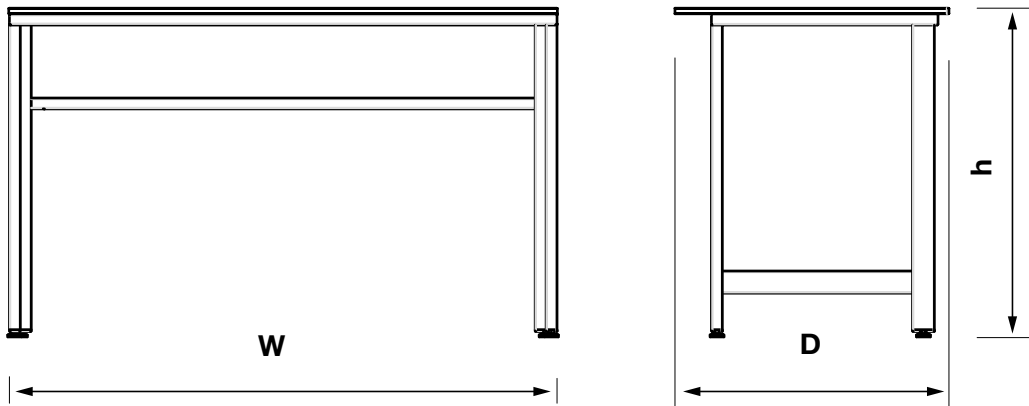

W (width, mm)	D (depth, mm)	h (worktop height, mm)
1.500	600	900
1.200	750	750
900	900	
750		

## MAX Load 200 Kg .....

Accessories and attachments such as electrical trunking brackets, side and back covers, or on demand new designs, can be attached at back beam and back legs nuts.


Primary Bench can be configured with appliances, check Primary Bench+ sheet for more info.

Check Primary Bench 5P sheet for info on extralarge sizes (up to 2.400 mm width)

Main frame fully welded steel, epoxy coated 9010 (check other options) with M5 nuts below worktop (back beam) and M6 at back legs for accesories and attachments.

Worktop material on demand. Advised HPL 16 mm high chemical resistance.

Feet, articulated with rubber base.

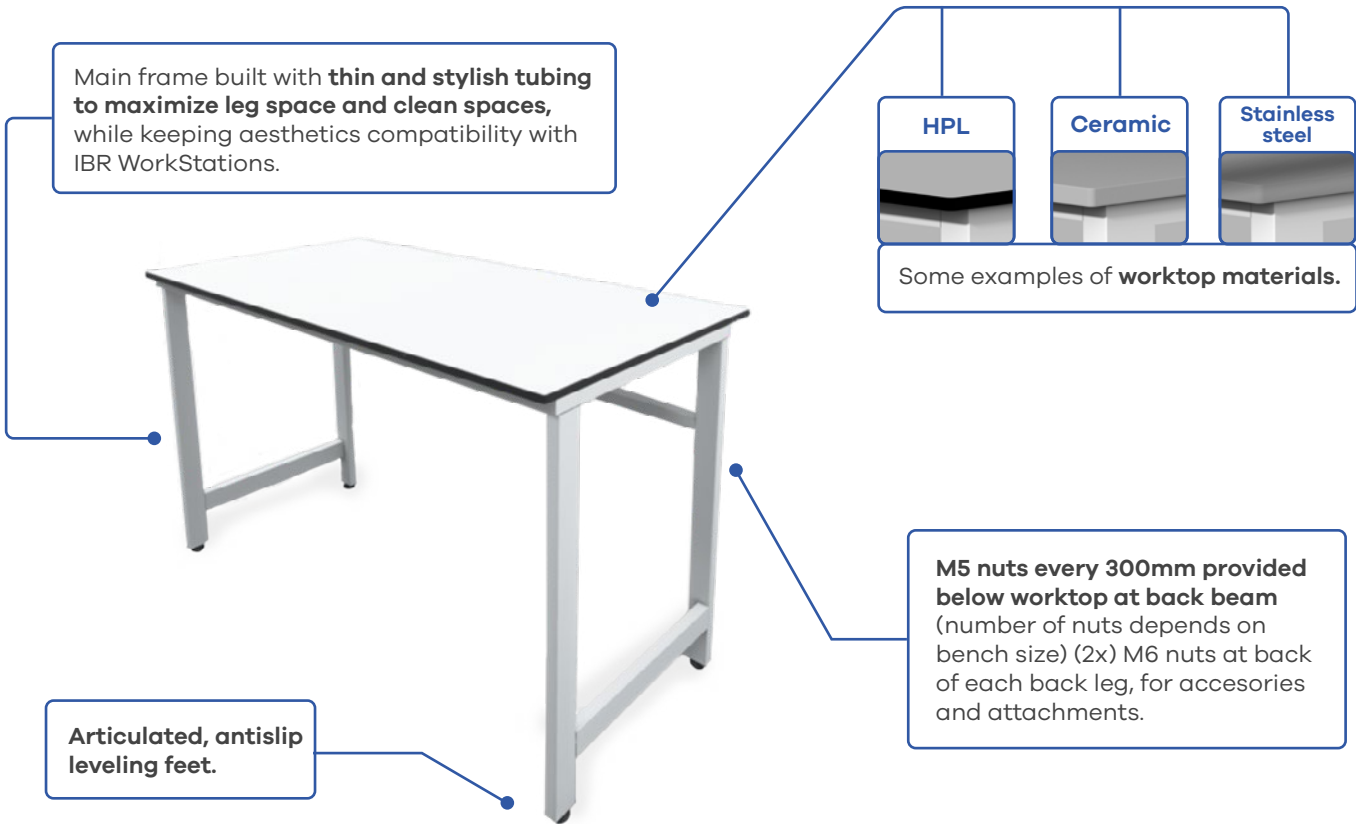

Primary Benches can also form islands and peninsulas on their own.



Primary Bench id an ideal aolution for a traditional approach in fixed positions, small rooms, corners, etc.



Makes a great complement for a comb layout setting the wall fixed benching as a starting point for comb's peninsulas, marking widths and spacing perfectly, by combining 1.500 mm width Primary Benching with 750 mm deep WorkStations (highlighted in green).

