

PRIMARY SINK

*Compact, practical
and ready for any configuration.*

.....




Primary Sink

Compact, practical and ready for any configuration.



Scan the QR code to see this product on our website or [click here](#)

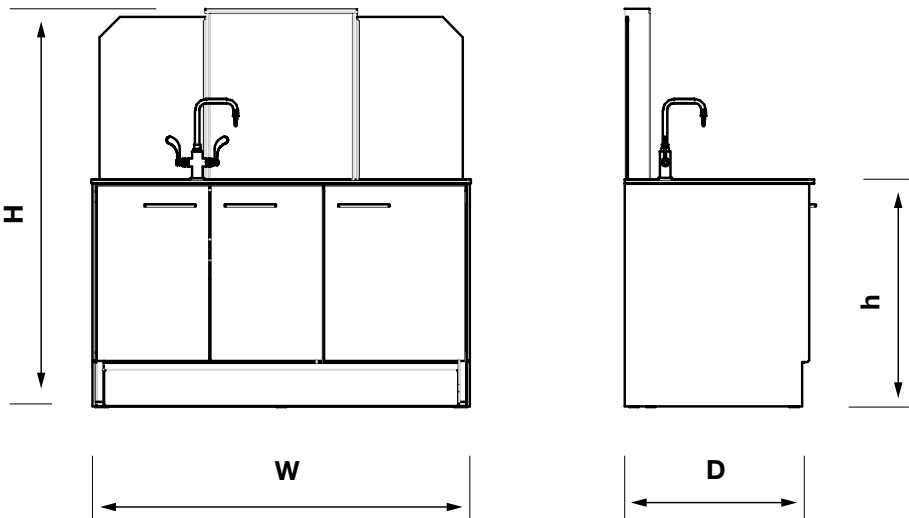

W (width, mm)	Doors (per width)	D (depth, mm)	h (worktop height, mm)	H (column height)
1500	3 (1.500)	600		
1200	2 (1.200)	750	900	1.570
900	2 (900)	900		UP TO CEILING
600	1 (600)			

.....
Fluids can be installed on worktop or at column front face (if present).


Designed for a quick installation on site.

A simple electrical setup can be installed at column on demand.

- | **Main core frame** fully welded steel, epoxy coated 9010 (check other options).
- | **Columns** (if present) manufactured in **metal sheet steel tube bolted to main frame, epoxy coated 9010** (Check other options) Size on demand.
- | **Outer panelling on demand** (HPL, Melamine, Metal sheet).
- | Worktop material **on demand**.
- | Sink material **on demand**.
- | **Feet, articulated with rubber base**.
- | **Soft closing hinges opening 165°**.

| **Fluids on demand** (water mixer, eye shower, RO Water).

| **Electrical setup up to 2 power sockets** can be installed at column (if present). Additionally, **circuit breaker setup on demand**.

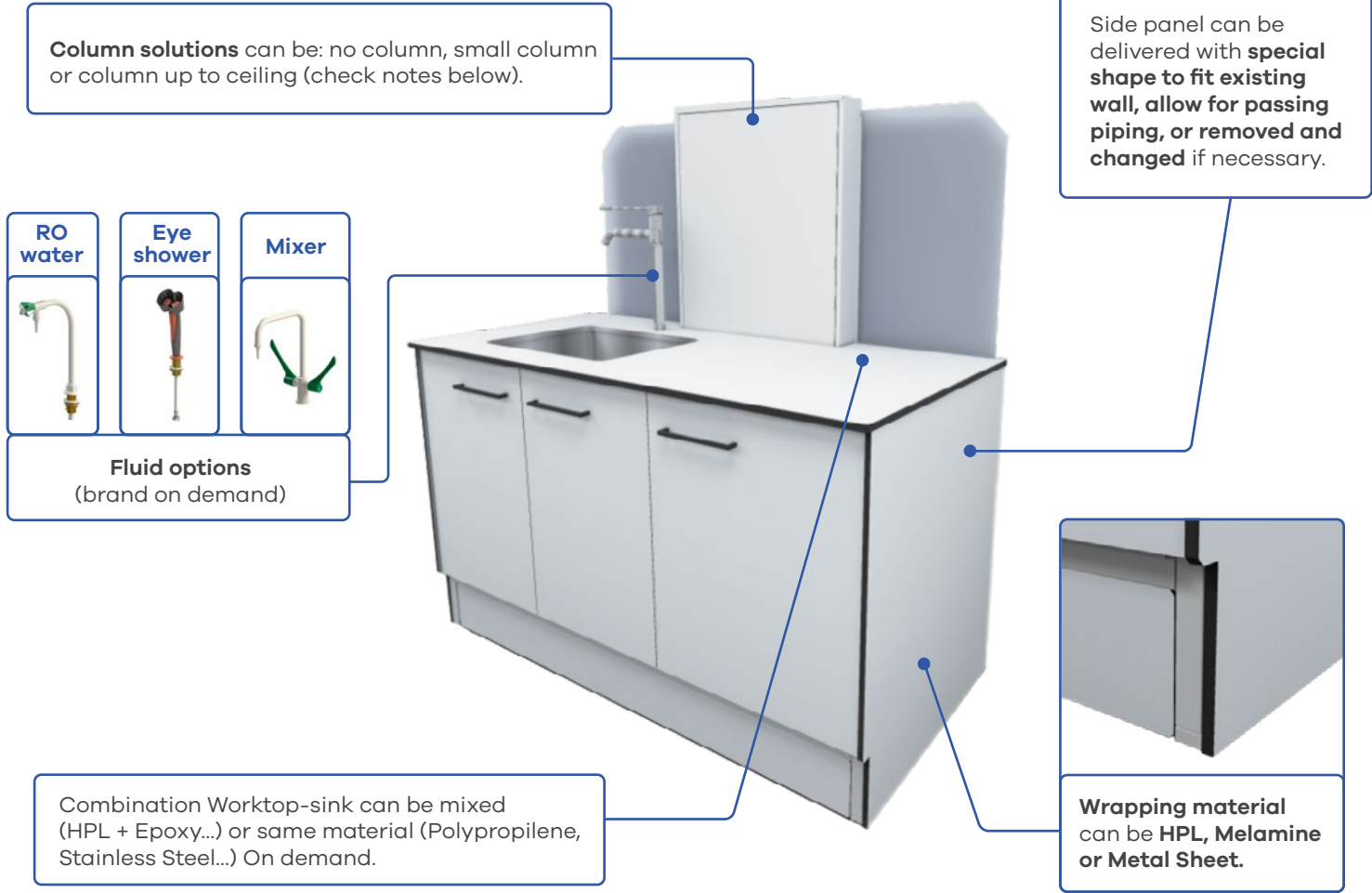

Clean setup for **concealed fluid hosing coming from ceiling**.



Same setup for island (covered back) or wall (naked back). Visible fluid hosing coming from the top.



Special setup with shelving installed and fluids coming from wall or floor.

Ask for special combinations if you need a unique setup

