

WORKSTATION BARCELONA

Our most iconic WorkStation;
Versatile. Distinctive. Refined.


.....



Our most iconic WorkStation;
Versatile. Distinctive. Refined.



Scan the QR code
to see this product
on our website or
[click here](#)

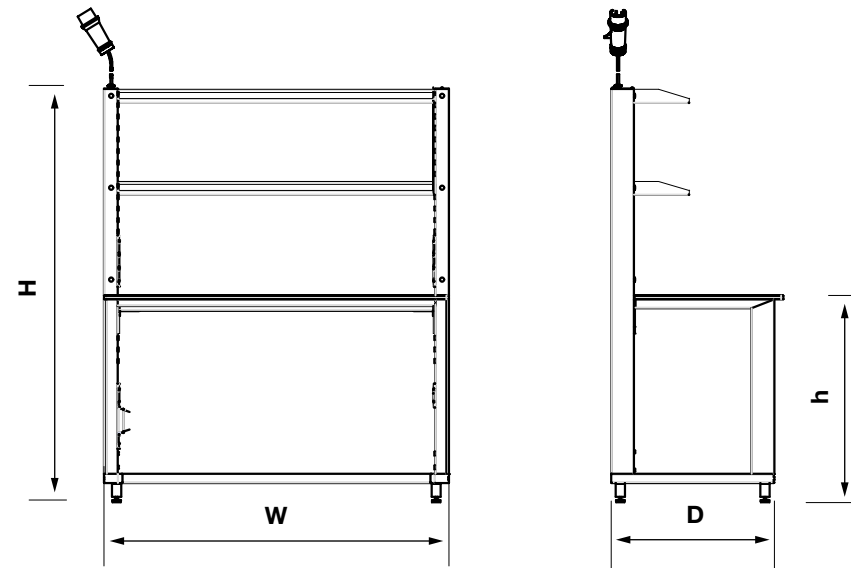

W (width)	D (depth)	h (worktop height)	H (column height)
1.500			
1.200	750	750	1.570
900	900	650	1.800
750			

MAX Load 350 Kg

Power supply side can be switched left to right and viceversa. Check **accessories pages** to see how shelves can be added, removed and changed position.

- **Underbench shelf** available.
- **Gas taps** can be attached to column.
- Some **other options** available.
- In addition, **special accessory attachments** can be design and manufactured on demand.

Main frame fully welded steel, epoxy coated 9010 (check other options) with M5 nuts below worktop for accesories and attachments.

Columns manufactured in steel tube bolted to main frame, epoxy coated 9010 (check other options), with **M8 nuts at internal faces** for accesories. **Nuts in front faces** as an option.

Worktop material on demand. Advised HPL 16 mm high chemical resistance.

Feet, articulated with rubber base, interchangeable with combi castor-foot.

Data setup on demand, advised to be magnet clip-routed.

Electrical setup 8 power sockets installed (4 each side above worktop). **Circuit breaker setup** on demand.

