

WORKSTATION OXFORD MNL


Simple, reliable & *precise*.

.....





Scan the QR code to see this product on our website or [click here](#)

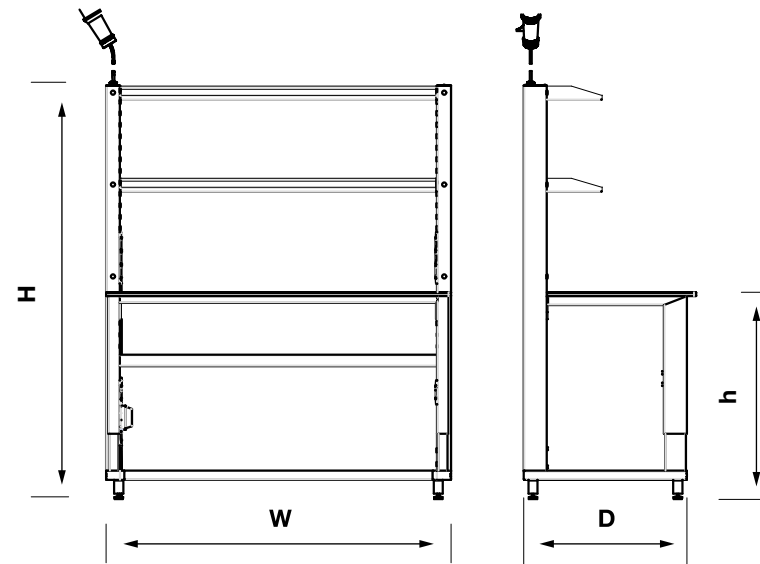

W (width, mm)	D (depth, mm)	h (worktop height, mm)	H (column height, mm)
1.500	750	SELECTABLE	1.570
1.200	900	750 - 800 - 850 - 900	1.800
900			

MAX Load 350 Kg

Power supply side can be switched left to right and viceversa. Check **accessories** pages to see how shelves can be added, removed and changed position.

- Underbench shelf available.
- Gas taps can be attached to column.
- Some other options available.
- In addition, special accessory attachments can be design and manufactured on demand.

Base frame fully welded steel, epoxy coated 9010 (check other options).

Worktop frame fully welded steel, to be bolted to base frame and columns at any of the selectable heights. **Provided with M5 nuts below worktop** for accesories and attachments.

Columns manufactured in steel tube bolted to main frame, **epoxy coated 9010** (check other options) with **M8 nuts at internal faces** for accesories. **Nuts in front faces as an option.**

Worktop material on demand. Advised HPL 16 mm high chemical resistance.







Feet, articulated with rubber base, interchangeable with combi castor-foot.

Data setup on demand, advised to be magnet clip-routed.


Electrical setup 8 power sockets installed (3 each side above worktop, 1 each side below worktop). **Circuit breaker** setup on demand .

Typical connection to ceiling options:

- CE connector plug
- UK plug
- Schuko plug
- Thorsman plug

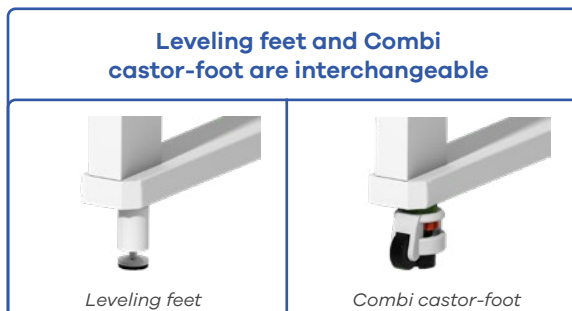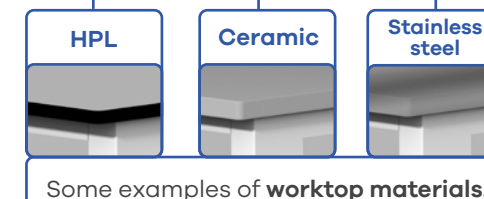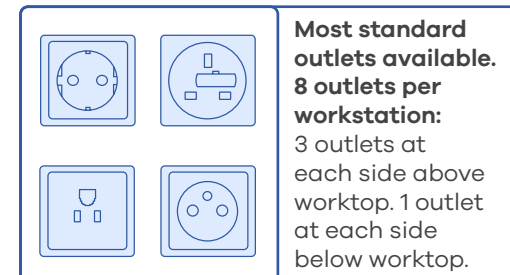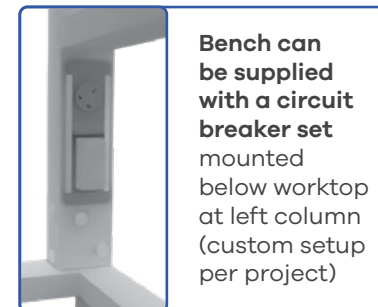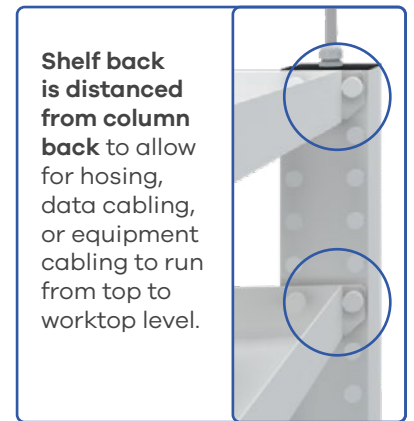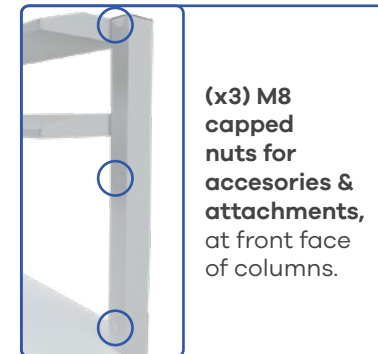

Worktop height can be changed by unbolting worktop frame from both telescopic front leg and column side, and bolting it again at any of the available positions.

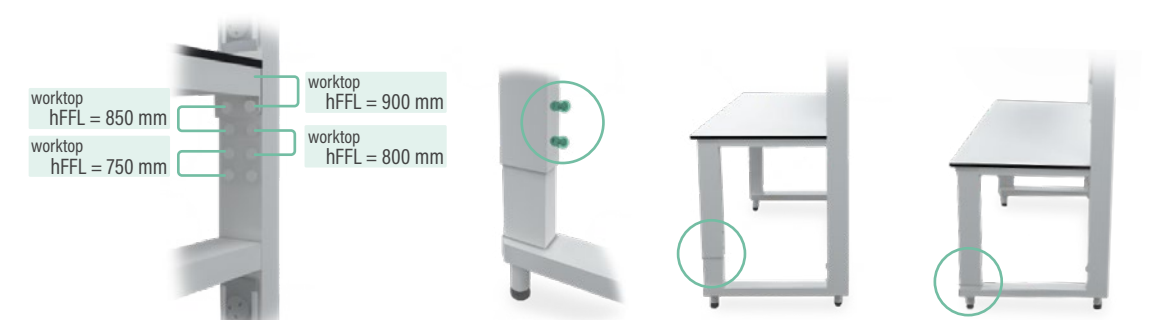

hFFL: high from floor level.

© IBR Innovación S.L. – All rights reserved. Partial or total reproduction of this catalog is strictly prohibited. See full legal notice in this document.