

PRIMARY + [T] 4P BENCH

All-in-one. Power, data and elevated service capacity. Ready for whatever comes next.

.....




Primary+ [T] 4P

All-in-one. Power, data and elevated service capacity. Ready for whatever comes next.



Scan the QR code to see this product on our website or [click here](#)

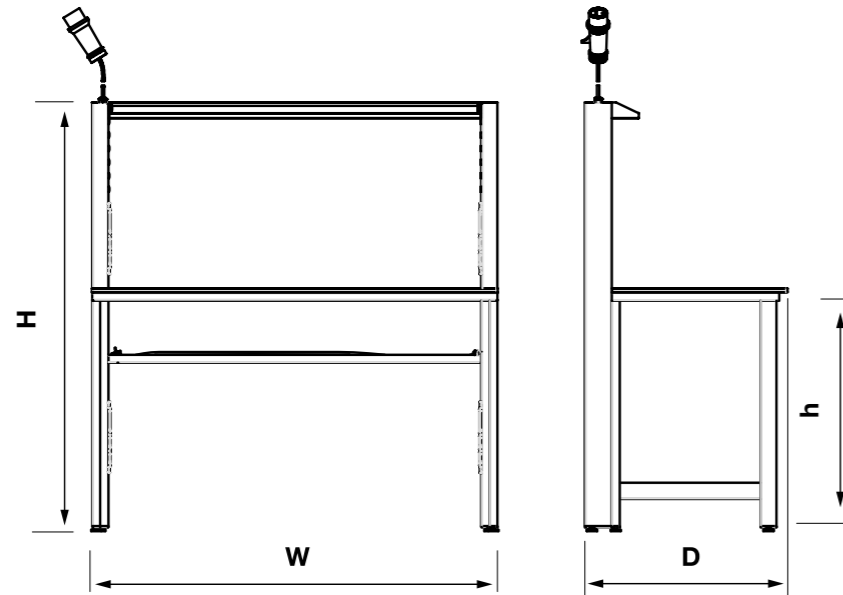

W (width, mm)	D (depth, mm)	h (worktop height, mm)	H (column height, mm)
1500	750 / 900	900 / 750	1600 / 1800
1200	750 / 900	900 / 750	1600 / 1800
900	750 / 900	900 / 750	1600 / 1800

MAX Load 350 Kg

Power supply side can be switched left to right and viceversa. Check **accessories pages** to see how shelves can be added, removed and changed position.

Primary Bench+ [T] 4P can always be built from a Primary Bench by adding the columns attachment and changing the worktop.

- **Gas taps** can be attached to column.
- Some **other options available**.
- In addition, **special accessory attachments** can be design and manufactured on demand.

Main frame fully welded steel, epoxy coated 9010 (check other options) with **M5 nuts below worktop** (back beam) and M6 at back legs for accesories and attachments.

Worktop material on demand. Advised HPL 16 mm high chemical resistance.


Feet, articulated with rubber base.

Bolted (detachable) set of electrical columns main body welded steel, epoxy coated 9010 (check other options).

Electrical setup: 8 power sockets installed (3 each side above worktop, 1 each side below worktop).

Circuit breaker setup on demand.

Data setup on demand, advised to be magnet clip-routed.

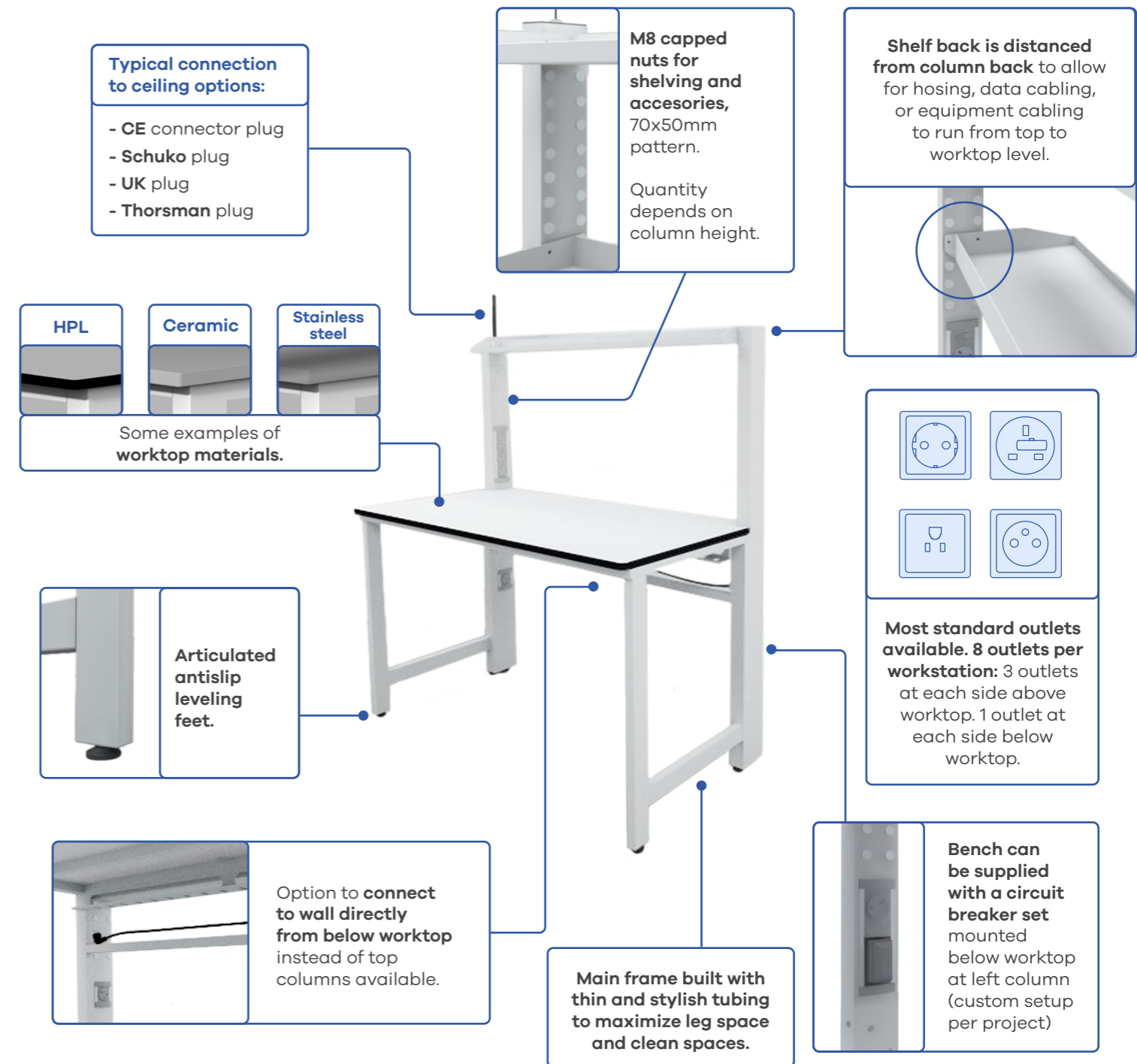

Primary+ [T] can work in combination with the other elements in the Primary System.

Also, **single column assembly** are available and fully working asymmetrical solutions.