

BALANCE BENCH

The timeless solution for precision and stability
—classic performance, discreet elegance.

.....




Balance

Minimalist and elegant—engineered for *rigidity, with a floating look.*



Scan the QR code to see this product on our website or [click here](#)

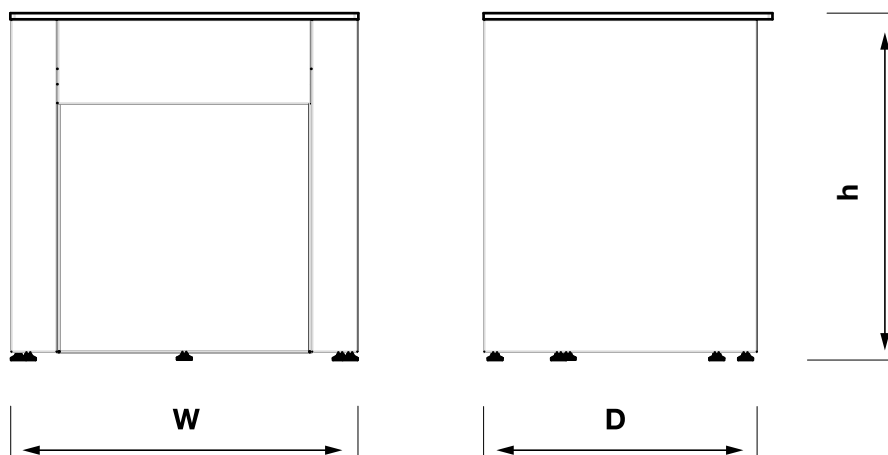

W (width, mm)	D (depth, mm)	h (worktop height, mm)
---------------	---------------	------------------------

900

750

900

Ask for special sizes on demand.



.....

Designed for quick installation on site, IBR Balance Bench is **delivered in 4 preassembled parts.** (5 if Leg Space Cover requested)

Custom sizes are available upon request.

Casing can be adapted from factory to fit special walls or skirting.

Internal frame fully welded steel, epoxy coated 9010 (check other options) **prepared with shock absorbant silent blocks at base,** and at top beams to receive weighed stoneware.


Weighed stoneware ~75kg provides **stability to core of balance bench,** sits on silent blocks at top of internal frame.

Casing manufactured in metal sheet, epoxy coated 9010 (check other options) **supplied preassembled with HPL worktop and x4 leveling feet** articulated with rubber base at its base.

Worktop material on demand. Advised HPL 16 mm high chemical resistance.

Leg space cover manufactured in metal sheet, epoxy coated 9010 (check other options) **supplied preassembled with HPL worktop and x4 leveling feet** articulated with rubber base at its base.

Weighing surface manufactured in black stoneware, smooth finished working surface.

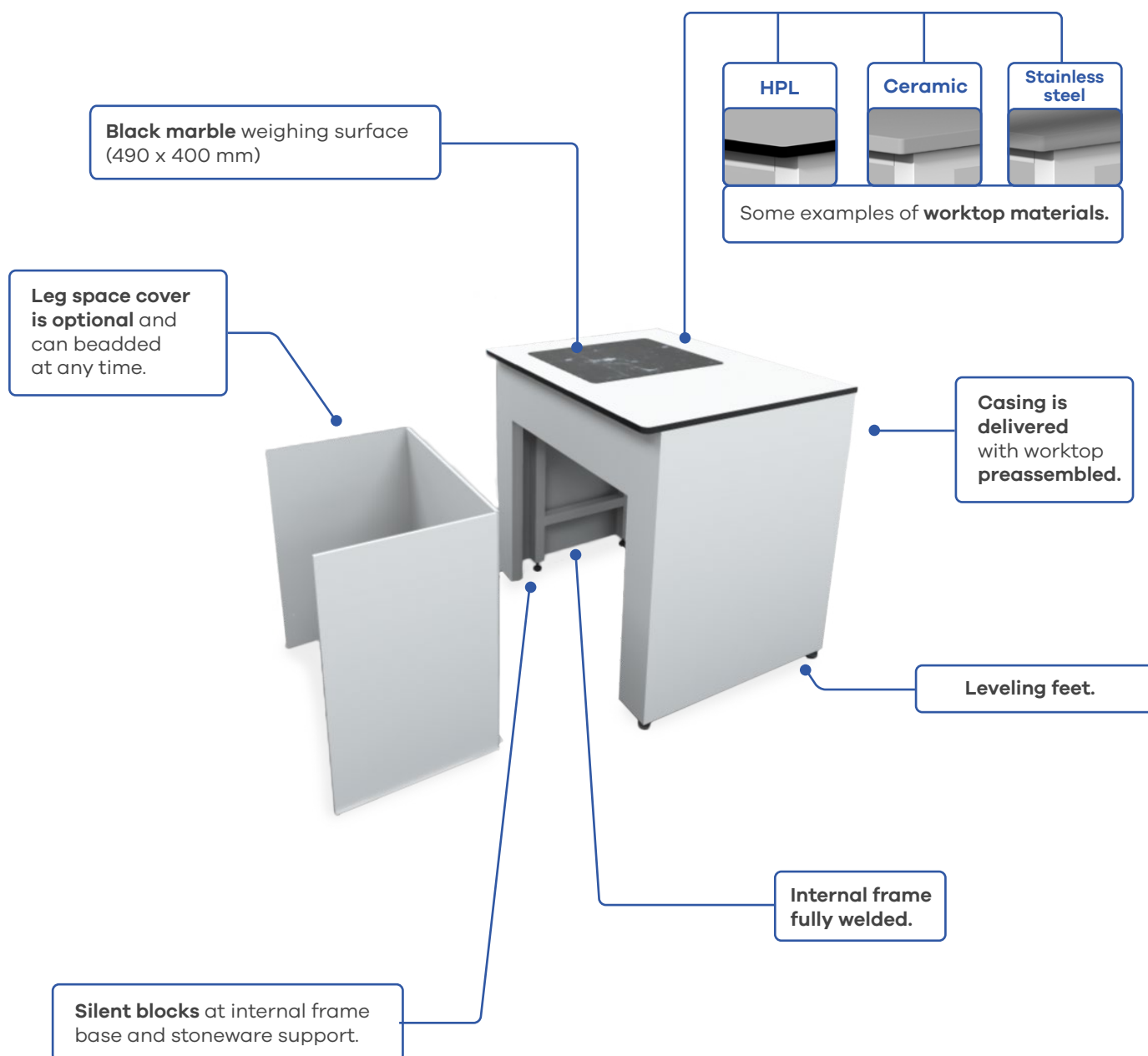

Technical Data Sheets

Visit our
website for
more info



www.ibr-innovasystem.com

Contact us at: hello@ibr-innovasystem.com

