

CANTILEVER WITH PANNELING

Minimalist and elegant—*engineered
for rigidity, with a floating look.*

.....




Cantilever with Panneling

Minimalist and elegant—*engineered for rigidity, with a floating look.*



Scan the QR code
to see this product
on our website or
[click here](#)

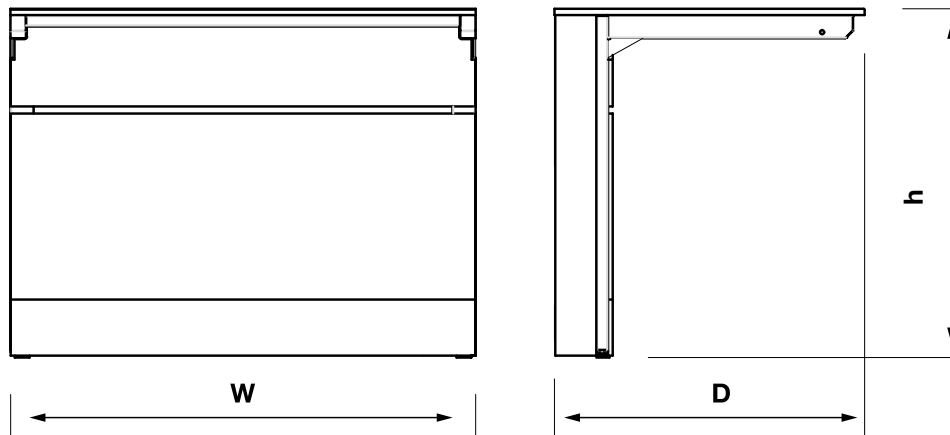

W (width, mm)	D (depth, mm)	h (worktop height, mm)
1200	600 / 750	900 / 750
900	600 / 750	900 / 750
600	600 / 750	900 / 750

Max recommended load - 200 kg

Check accesories pages to see **underbench attachments options** available.

In addition, **special accesory attachments** can be design and manufactured on demand.

Panelling can be manufactured on HPL, melamine or metal.


Main frame fully welded steel, epoxy coated 9010 (check other options) **with M5 nuts below worktop** for accesories and attachments.

Worktop material on demand. Advised HPL 16 mm high chemical resistance.


Panelling on demand (HPL, Melamine, Metal sheet). **Fixed panelling bolted to main structure.** **Accesible panelling magnet** attached to main structure.

Feet, articulated with rubber base.

Wall bracket should be attached to wall with (x3) plugs min. 200 kg pulling force each, or similar means depending on wall material.

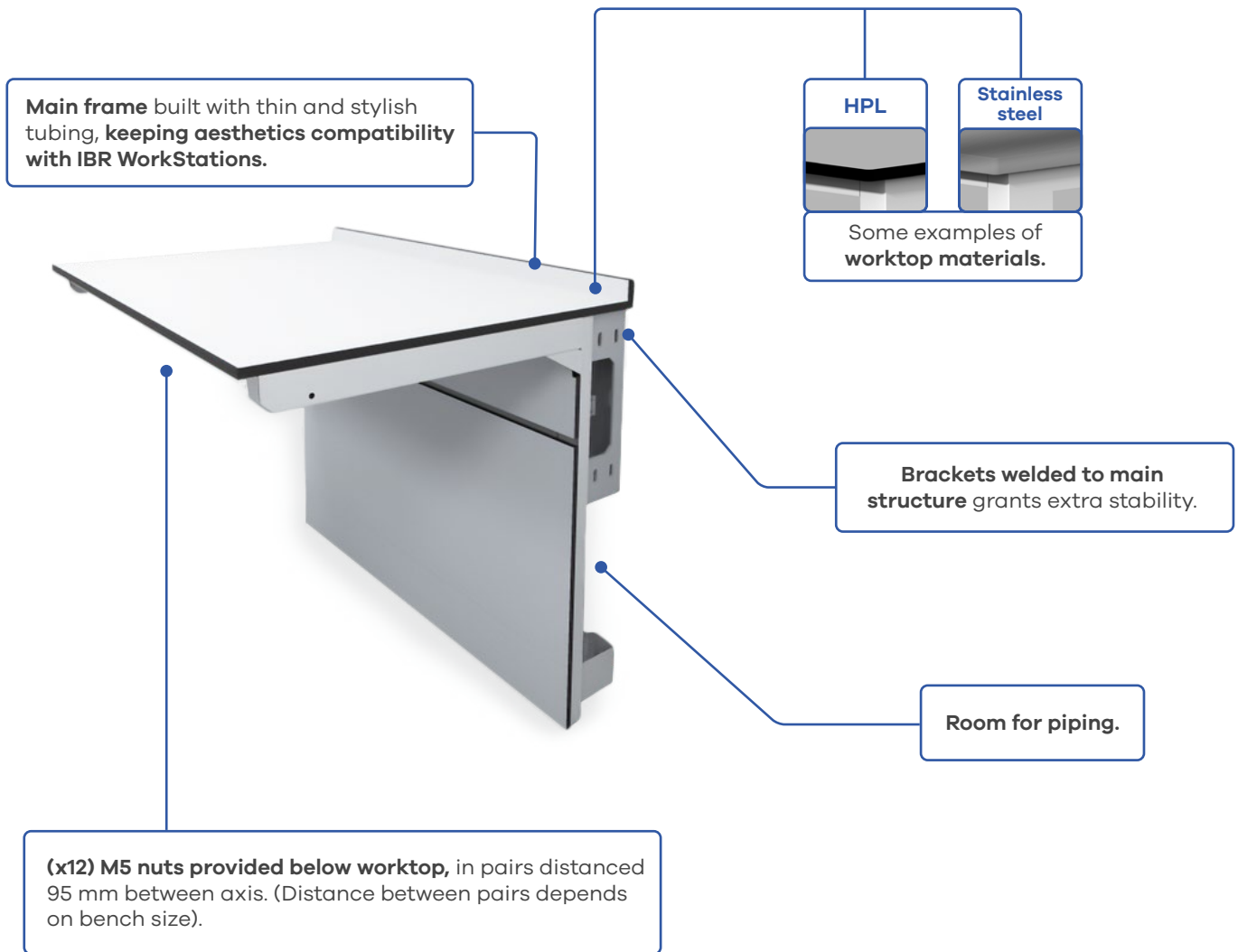

Side cover at the end of a row of Cantilever benches can be **shaped on demand** to fit wall skirting or piping. (Highlighted in green).

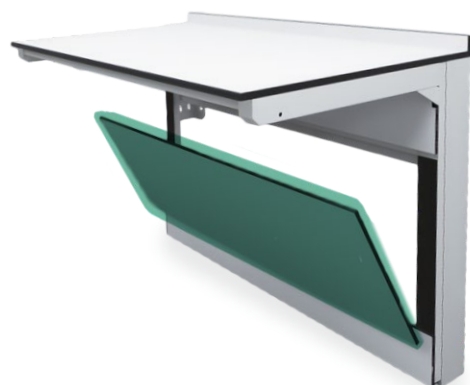

Front panel is attachable to main structure with magnets to grant easy accesibility. (Highlighted in green).

Technical Data Sheets

Visit our
website for
more info



www.ibr-innovasystem.com

Contact us at: hello@ibr-innovasystem.com

