

WORKSTATION OXFORD MNL


Simple, reliable & *precise*.

.....





Scan the QR code
to see this product
on our website or
[click here](#)

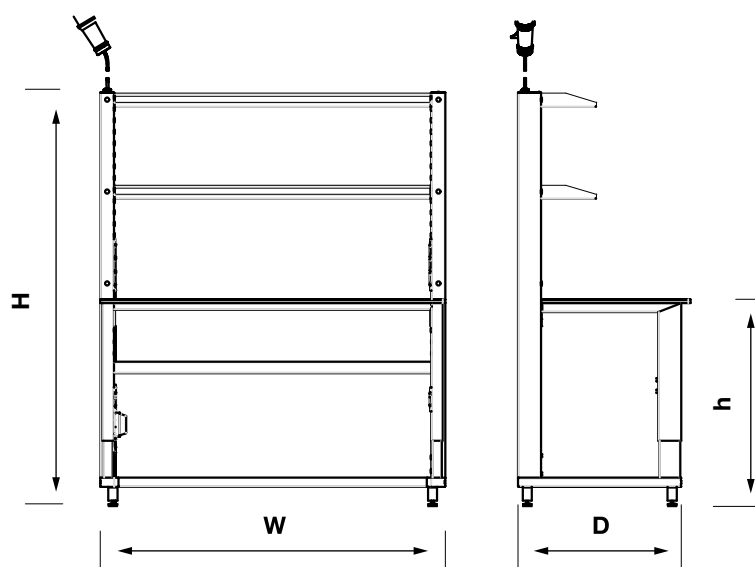

W (width, mm)	D (depth, mm)	h (worktop height, mm)	H (column height, mm)
1.500	750 / 900	SELECTABLE	1.570 / 1.800
1.200	750 / 900	750 - 800 - 850 - 900	1.570 / 1.800
900	750 / 900		1.570 / 1.800

Max recommended load - 400 kg

Power supply side can be switched left to right and viceversa. Check **accessories pages** to see how shelves can be added, removed and changed position.

- **Underbench shelf** available.
- **Gas taps** can be attached to column.
- Some **other options** available.
- In addition, **special accessory attachments** can be design and manufactured on demand.

Base frame fully welded steel, epoxy coated 9010 (check other options).

Worktop frame fully welded steel, to be bolted to base frame and columns at any of the selectable heights. **Provided with M5 nuts below worktop** for accesories and attachments.


Columns manufactured in steel tube bolted to main frame, **epoxy coated 9010** (check other options) with **M8 nuts at internal faces** for accesories. **Nuts in front faces as an option.**

Worktop material on demand. Advised HPL 16 mm high chemical resistance.

Feet, articulated with rubber base, interchangeable with combi castor-foot.


Data setup on demand, advised to be magnet clip-routed.

Electrical setup 8 power sockets installed (3 each side above worktop, 1 each side below worktop). **Circuit breaker** setup on demand .

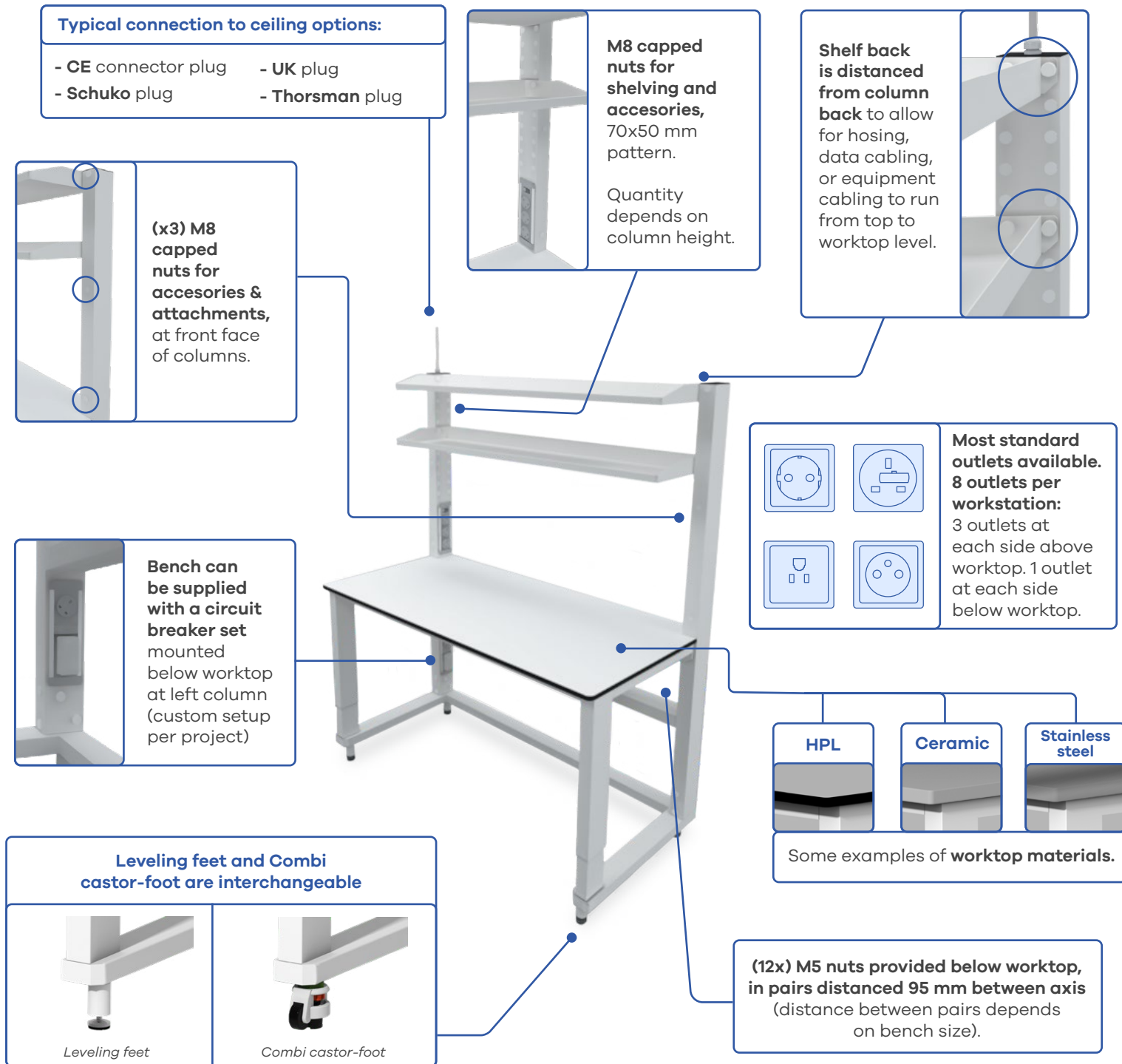

Worktop height can be changed

by unbolting worktop frame from both telescopic front leg and column side, and bolting it again at any of the available positions.

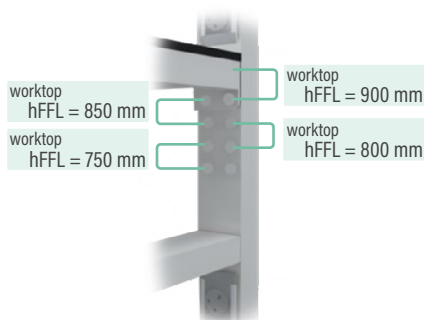

hFFL: high from floor level.

© IBR Innovación S.L. – All rights reserved. Partial or total reproduction of this catalog is strictly prohibited. See full legal notice in this document.

Technical Data Sheets

Visit our
website for
more info



www.ibr-innovasystem.com

Contact us at: hello@ibr-innovasystem.com

